




Medial

Innovating Healthcare, Empowering Life

Disposable Products Catalogue

www.medial.com.tr

disposet

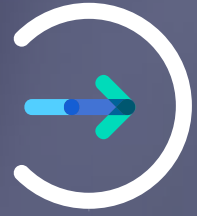

PRESSURE MONITORING SET



We offer a comprehensive range of Pressure Monitoring Set designed to meet the diverse needs of healthcare professionals. Our product lineup includes both wired and wireless options, providing flexibility to suit different clinical environments and preferences. Our sets are available in two key configurations as Invasive Sets and Safe Set.


PRESSURE MONITORING SET

INVASIVE-STANDARD

Standard Pressure Monitoring Set

Pressure monitoring set converts blood pressure into electrical signals via transducers, and these signals are transmitted to a monitor either using a direct cable from the transducer or a cable from the wireless contact board. Pre-calibrated sensors with constant electrical parameters enable reliable and continuous follow-up of the patient's blood pressure. In the wired models, a jack extends from the transducer body and a cable is then connected directly to the monitor. In the wireless models, the transducers are placed on a contact board, and data is transmitted to the monitor by a cable connecting the monitor to this board. The contact board is reusable and requires a single cable.

- Easy flush
- Available with single, double, triple and quadruple transducers
- Wireless data transfer (optional)
- Dielectric gel barrier

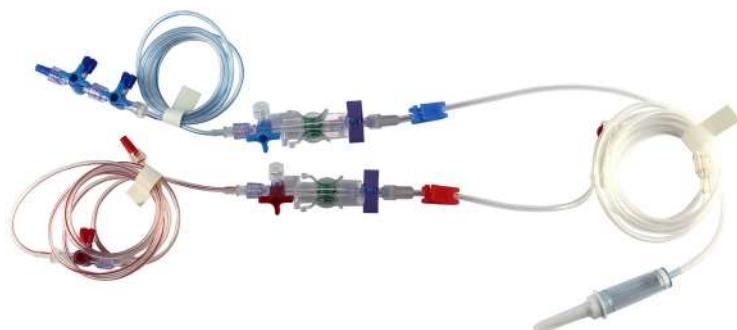


SAFE SET

Safe Set Pressure Monitoring Set

Safe Set blood pressure monitoring set, incorporating a closed blood sampling system, provide all the advantages of conventional pressure monitoring system while also featuring an advanced blood sample collection mechanism. This innovative system enables the collection of blood samples without the need for a needle and ensures protection from atmospheric exposure. Consequently, the risk of contamination and needle-related injuries is significantly reduced.

- Reliable and continuous blood pressure monitoring
- Easy flush
- Wireless data transfer (optional)

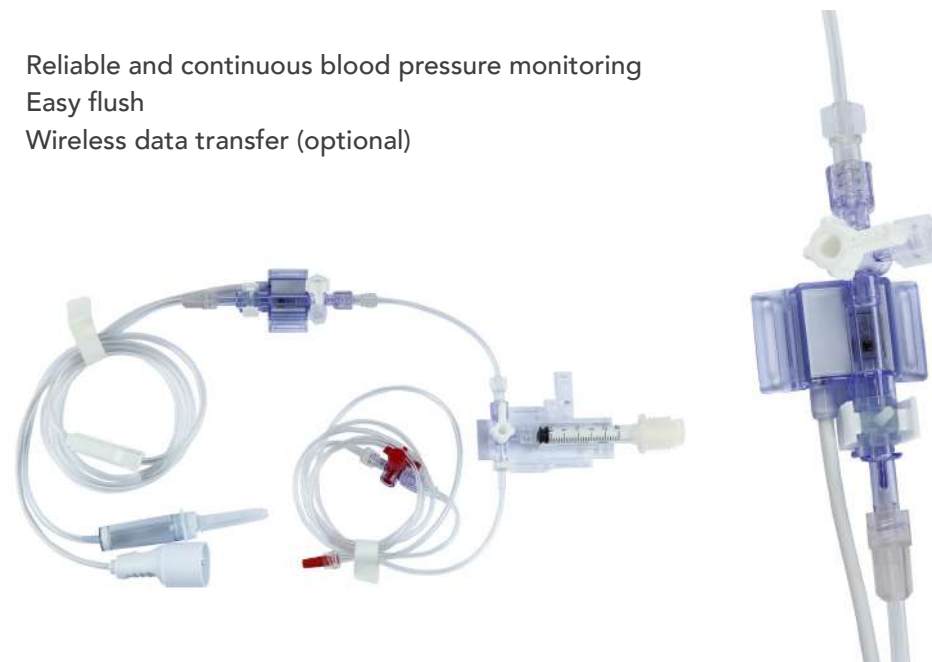

WIRED TRANSDUCER

Transducer Body – With Cable Connection

The transducer is designed to convert blood pressure into electrical signals and features a cable connection port integrated directly into the transducer body. This port connects to a monitor via a cable, and these connections are available in various configurations. Transducer bodies can be provided without a flush unit, if requested.



WIRELESS TRANSDUCER

Transducer Body – The transducer, which converts blood pressure into electrical signals, is linked to a dedicated board that transmits these signals wirelessly. The transducer body is designed without a direct cable connection, using a dedicated board to wirelessly transmit signals to the monitor via a connecting cable for display.



CLOSED BLOOD SAMPLING

Closed blood sampling components are essential for pressure monitoring set. They enable the secure collection of blood samples while preserving the integrity of the set and reducing risks associated with environmental contamination.





Innovating Healthcare Empowering Life

Address:

Bestekar Şakir Ağa Sok. Trump Villa, 6/1 Şişli - İstanbul / Turkey

E-Mail:

export@medial.com.tr

www.medial.com.tr